

## Spare Parts

### PRESTO AQUA LUX

MODEL	PRODUCTS	ORDER CODE
PRESTO AQUA LUX Control Unit		Z1038A002
PRA-LUX-HP PRESTO AQUA LUX Handpiece		T1002
PR-AQL Hose Coupling Hose		Z1050

MODEL	PRODUCTS	ORDER CODE
PR-AQL03 PRESTO AQUA LUX Cartridge		T1002250
Handpiece Stand		Z095206
AC Adaptor PR-AQL		U1044

MODEL	PRODUCTS	ORDER CODE
PRESTO AQUA LUX Cord (120V)		40010301
PRESTO AQUA LUX Cord (230V)		U439550
PTL O-ring Set (Pack of 5)		Y900580

### PRESTO AQUA II

MODEL	PRODUCTS	ORDER CODE
PRESTO AQUA II Unit		Z308015
PRA-HP PRESTO AQUA Handpiece		T816

MODEL	PRODUCTS	ORDER CODE
Hose for PRESTO AQUA II		Z308050
QD-J B3 Coupling		T354051

MODEL	PRODUCTS	ORDER CODE
PR-AQ03 PRESTO AQUA Cartridge		T816016

### PRESTO AQUA LUX / PRESTO AQUA II

MODEL	PRODUCTS	ORDER CODE
AFC-40 Foot Control		Z127
Chuck		T815019
Secondary Hose for PRESTO (2 m)		Z101240
Water Filter Set (water tube x 2, water filter, water filter case)		U387040
Water Filter with Case		Y900069

MODEL	PRODUCTS	ORDER CODE
Water Filter		U387042
Tube Ring (Includes 2)		Y900077
PRA II Bottle Set		Z308040
PRA II Bottle Lid		Z308041
PRA II Bottle		Z308601

MODEL	PRODUCTS	ORDER CODE
PRA II Bottle O-ring set (Pack of 2)		Y900327
Bracket Set		Z308060
Chuck Wrench		C034751
Bur-inserting Tool		T815030

### PRESTO II

MODEL	PRODUCTS	ORDER CODE
PRESTO II Unit		Z307015
PR-HP PRESTO Handpiece (Includes 2m hose)		T815
AFC-40 Foot Control		Z127
Handpiece Stand		Z010134

MODEL	PRODUCTS	ORDER CODE
PR-03 PRESTO Cartridge		T815016
Chuck		T815019
Hose for PRESTO (2 m)		T815220
Secondary Hose for PRESTO (2 m)		Z101240

MODEL	PRODUCTS	ORDER CODE
LAB Spanner Wrench		H005751
Chuck Wrench		C034751
Bur-inserting Tool		T815030

# NSK

## LUBE FREE LAB AIR TURBINES

PRESTO AQUA LUX  
PRESTO AQUA II  
PRESTO II

# optimum performance & precision

## LED

Creating the Perfect Fit with LED Illumination

<b>NAKANISHI INC.</b> www.nsk-inc.com 700 Shimohinata Kanuma-shi Tochigi 322-8666, Japan TEL : +81 (0)289-64-3380 FAX : +81 (0)289-62-5636	<b>NSK Europe GmbH</b> www.nsk-europe.de Elly-Beinhorn-Strasse 8 65760 Eschborn, Germany TEL : +49 (0)6196 77606-0 FAX : +49 (0)6196 77606-29	<b>NSK United Kingdom Ltd</b> www.nsk-uk.com Office 5, Gateway1000, Arlington Business Park, Whittle Way, Stevenage, SG1 2FP, UK TEL : +44 (0)1438 310670 FAX : +44 (0)1438 310671	<b>NSK Oceania Pty Ltd</b> www.nsk-inc.com Unit 22, 198-222 Young St. Waterloo, Sydney, NSW 2017, Australia TEL : +61 (0)2 8306 3000 FAX : +61 (0)2 8306 3030	<b>NSK America Corp</b> www.nskamerica.com 700 Cooper Court Schaumburg, IL 60173, USA TEL : +1 888 675 1675 FAX : +1 800 838 9328	<b>NSK Middle East</b> www.nsk-inc.com Room 6EA-701, 7th Floor, East Wing No.6 Dubai Airport Free Zone, PO Box 54316 Dubai, UAE TEL : +971 (0)4-204-5507 FAX : +971 (0)4-204-5508	<b>NSK Asia</b> www.nsk-inc.com 1 Maritime Square, #09-33 HarbourFront Centre, Singapore 099253 TEL : +65-6376-5018 FAX : +65-6376-0233



# PRESTO AQUA LUX

Lube Free Air Turbine Handpiece with LED and Water Spray

PRESTO AQUA LUX enhanced functionality with LED light has been introduced as part of the new line up of the PRESTO AQUA product range used for a wide range of dental laboratory applications: from precision cutting of high-strength ceramics to working with alumina and zirconia, requiring a water spray coolant.


The PRESTO AQUA LUX with LED is a single device that covers all of the delicate and varied operations in the dental laboratory, assuring high performance and adjustability.

## Advantages of using LED

**Clearer:** LEDs generate natural daylight quality light which assures clearer vision when compared to halogen light.

**Durable:** LED light sources are longer lasting and more durable than halogen bulb.

**Safer:** LED light is safer as it generates less heat than a halogen bulb. Natural daylight quality LED light is easy on the eye.

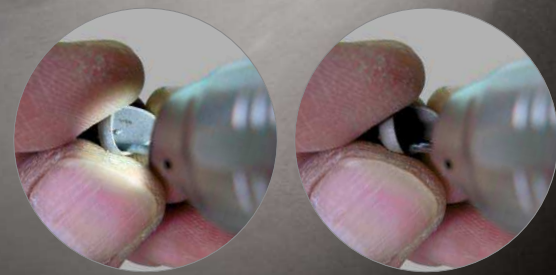

LED light

Non-Optic

### LED light enables you to work more efficiently

- Accurate lighting at the point of operation does not hide any detail, allowing you to focus on your work more easily.
- NSK LED is equivalent to daylight: it is easy on the eye, reflecting true coloring.
- LED illuminates fine unevenness and air bubbles of the internal frames which can not be illuminated under the normal bench lighting, thus dramatically improving the quality of internal adjustments.
- Enables fine recontouring effectively.
- Enables precise work and materials assessment under the microscope, which improves efficiency and reduces fatigue.

## Easy adjustment of the water coolant spray, LED light and air spray

It is possible to adjust the intensity of the LED, water coolant flow and spray air seamlessly. This adjustability is particularly useful when working with a wide range of materials and allows you to choose the most appropriate setting for your work at any time.




### ● Adjustment of LED light

It is possible to adjust the LED light in order to prevent light reflection and to set the most appropriate light volume for particular working conditions.

### ● Water Coolant Options

The water coolant can be seamlessly adjusted from water droplets to a fine water/air mist spray. Both the volume of water and the spray air can be adjusted independently to achieve the ideal working conditions for any material and application: very fine minimal water spray used for preventing micro-cracks to higher volume of water spray used for flushing away debris during cutting.

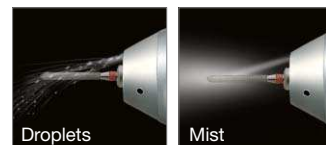

# PRESTO AQUA LUX **Optic**

Lube Free Air Turbine Handpiece with LED and Water Spray



MODEL PR-AQL Hose ORDER CODE Z1050  
• LED incorporated

MODEL PRA-LUX-HP ORDER CODE T1002

- 320,000 min<sup>-1</sup>
- Coolant solution spray
- Requires no lubrication
- Swivel twist-free handpiece
- Low noise and vibration
- Unique dust-proof mechanism

MODEL PRESTO AQUA LUX (120V) ORDER CODE Y1001199

MODEL PRESTO AQUA LUX (230V) ORDER CODE Y1001151

- Contents • PRESTO AQUA LUX Control unit • PRESTO AQUA LUX handpiece • PR-AQL Hose  
• Foot control (AFC-40) • Secondary hose for PRESTO • Handpiece stand  
• PRESTO AQUA LUX Cord (120V/230V) • AC Adaptor

- Simple and easy removal of the water bottle
- Seamless adjustments of LED light intensity, water volume and spray air
- On-Off function of LED light and water coolant spray • Capable of running independently with the incorporated water bottle or connected to mains water supply

# PRESTO AQUA II **Non-Optic**

Lube Free Air Turbine Handpiece with Water Spray



MODEL QD-J B3 ORDER CODE T354051  
• PRESTO AQUA coupling

MODEL PRA-HP ORDER CODE T816

- 320,000 min<sup>-1</sup>
- Coolant solution spray
- Requires no lubrication
- Swivel twist-free handpiece
- Low noise and vibration
- Unique dust-proof mechanism

MODEL PRESTO AQUA II ORDER CODE Y150023

- Contents • PRESTO AQUA II unit • PRESTO AQUA handpiece • QD-J B3 coupling  
• Hose for PRESTO AQUA II • Foot control (AFC-40) • Secondary hose for PRESTO

- Simple and easy removal of the water bottle
- Seamless adjustments of water volume and spray air • On-off function of water coolant spray
- Capable of running independently with the incorporated water bottle or connected to mains water supply

# PRESTO II **Non-Optic**

Lube Free Air Turbine Handpiece



MODEL PR-HP ORDER CODE T815

- 320,000 min<sup>-1</sup>
- Requires no lubrication
- Low noise and vibration
- Unique dust-proof mechanism
- Lightweight and compact

MODEL PRESTO II ORDER CODE Y150022

- Contents • PRESTO II unit • PRESTO handpiece • Handpiece stand  
• Hose for PRESTO • Foot control (AFC-40) • Secondary hose for PRESTO